

CDR Systems
Precision Patient Positioning Products

www.cdrrsys.ca info@cdrrsys.ca
403.850.7035 1.855.856.7035

FREEDOM TOTAL BODY PRECISION PATIENT IMMOBILIZATION SYSTEM, FOR SRS, SBRT AND MORE


PRODUCT OVERVIEW

A complete SRS/SBRT immobilization platform for both Cranial and Extracranial Stereotactic immobilization with command of pitch and roll axis of adjustment for the head, neck and shoulders

The Freedom System is a total solution developed with clinical input from radiation therapy centers world wide

Arms up Configuration shown in picture

ARMS UP CONFIGURATION

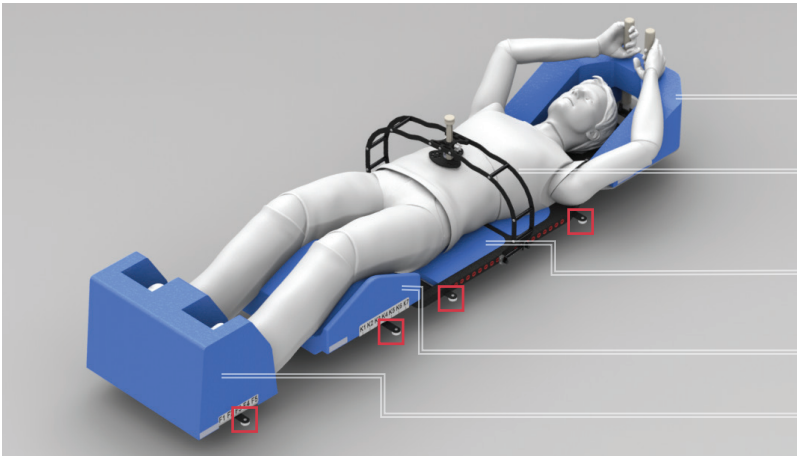


- Carbon fiber indexing bar set (match to couch top)
- ▶ Indexable Module for Arms Up/Back position
- ▶ Indexable Compression Bridge
- ▶ Filler Modules | Integrated table pads, removable for Indexable Vacuum Lock Cushion
- ▶ Indexable Module for Leg & Knee Positioning
- ▶ Indexable Module for Foot Positioning

Common setup uses:

SBRT patient immobilization with abdominal compression. Conventional lung. Supine breast. Pelvic treatment best suited to arms up position

Engineered for extremely low density, with a weight of only 1.7kg (3.8lbs) and easy to storage

Compatible with all treatment couches including: Varian, Elekta, Siemens, Tomotherapy and others.

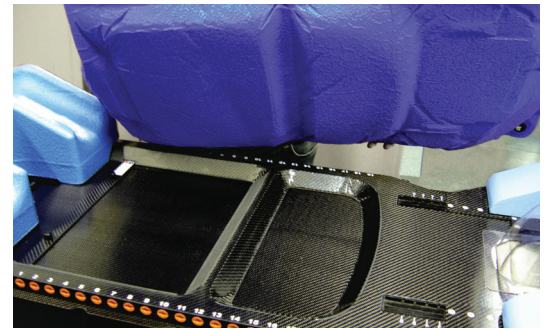


Vacuum cushion indexed automatically by the specifically engineered contours of the Freedom system



Freedom System Arms Up Configuration, shown with Vacuum Lock Bag

ARMS DOWN CONFIGURATION

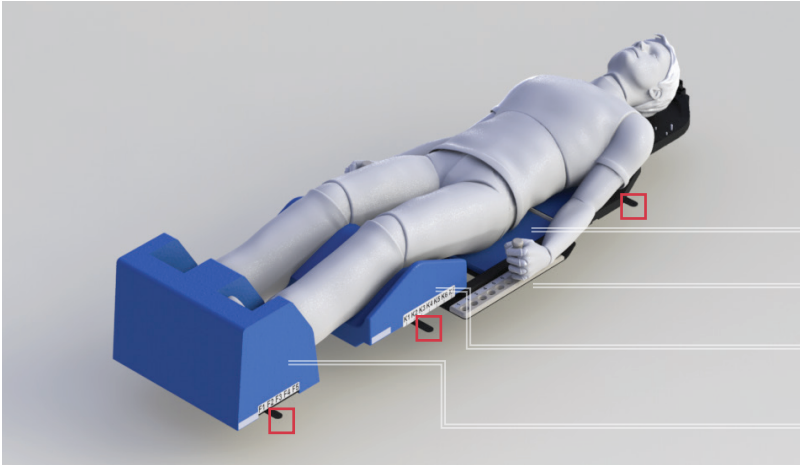

- Carbon fiber indexing bar set (match to couch top)
- ▶ Filler Modules | Integrated table pads, removable for Indexable Vacuum Lock Cushion
- ▶ Indexable Handles
- ▶ Indexable Module for Leg & Knee Positioning
- ▶ Indexable Module for Foot Positioning

Common setup uses:
 Frameless SRS for spherical tumors not requiring
 ** pitch & roll adjustability IGRT, IMRT head, neck
 and shoulder immobilization. Supine Prostate


**Optional (ID# PH08) for prone head and neck setup
 **For pitch and roll see CDR couch extension options*

*Engineered for extremely low density, with a weight of only 1.7kg (3.8lbs)
 and easy to storage*

*Compatible with all treatment couches including: Varian, Elekta,
 Seimens, Tomotherapy and others.*




Prone Head and Neck Immobilization System (Code: PH08)

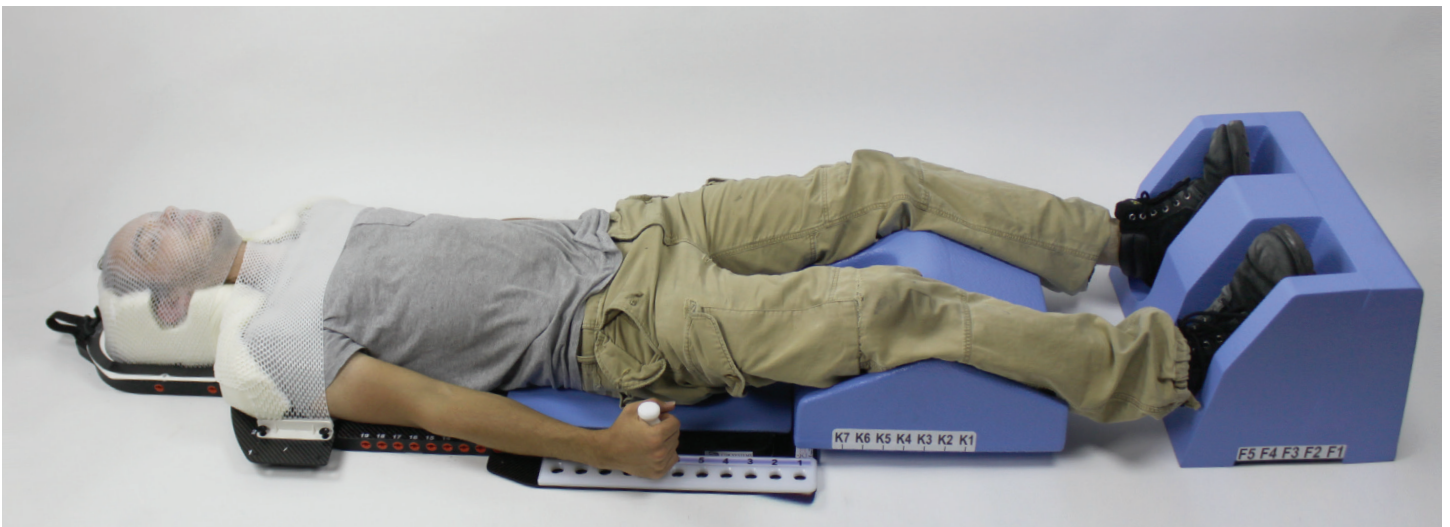

Freedom System Arms Down Configuration, shown with Intuition Mold, Shoulder Tech and Full LT Thermoplastic Mask

PARTS AND ACCESSORIES

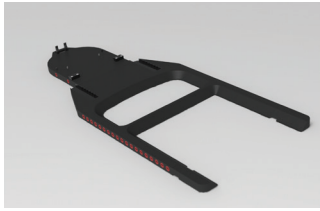

LB-SBRT-1

Freedom system indexable overlay


MR For MR safe version add "-K" to the order code



STA-12

ShoulderTech™ module compatible with: LB-SBRT, SB02 and LB02

MR For MR safe version add "-K" to the order code

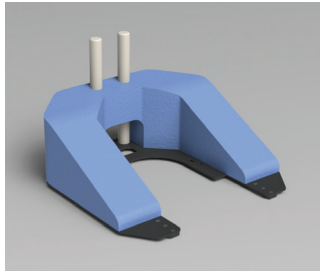


AU-3

Indexable module for Arms up configuration

AU-3-SA

Stand alone indexes directly to couch top


MR For MR safe version add "-K" to the order code

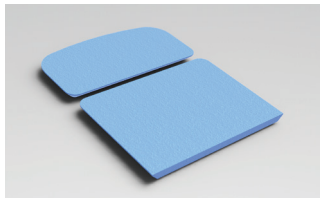

SBRT-FILL

Integrated table pad/filler modules

MR

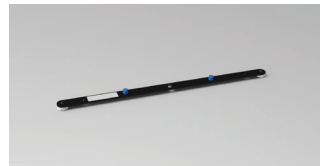

IB06-CF

Carbon fiber radiotranslucent indexing bar set of 4

MR For MR safe version add "-K" to the order code



SBRT-B1

Large indexable carbon fiber abdominal compression Bridge (shown with large compression post assembly CMP-L)


MR For MR safe version add "-K" to the order code



SBRT-B2

Small indexable carbon fiber abdominal compression Bridge (shown with small compression post assembly CMP-S)

MR For MR safe version add "-K" to the order code

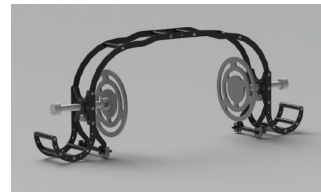

SBRT-B2-SP1-AC

Small indexable carbon fiber Bridge configured with index side paddles and arms down support for Arms Down patient positioning (available on large: SBRT-B1-SP1-AC)


MR For MR safe version add "-K" to the order code



SBRT-B1-SP1

Large indexable carbon fiber Bridge configured with side paddles (available on small bridge SBRT-B2-SP1)


MR For MR safe version add "-K" to the order code

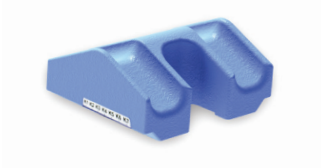

KN1-WDG-1

Indexable module for Knee and Leg positioner


MR

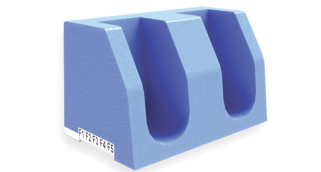

F-WDG-I

Indexable module for Foot and Leg positioner

MR

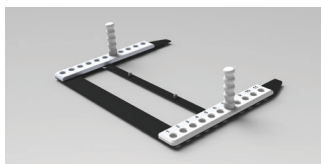

LB-HP

Arms down hand indexing module

MR For MR safe version add "-K" to the order code



LB-SBRT-WM

LB-SBRT couch overlay storage/ wall mount

MR For MR safe version add "-K" to the order code