




ADVANCED PRONE BREAST PATIENT POSITIONING

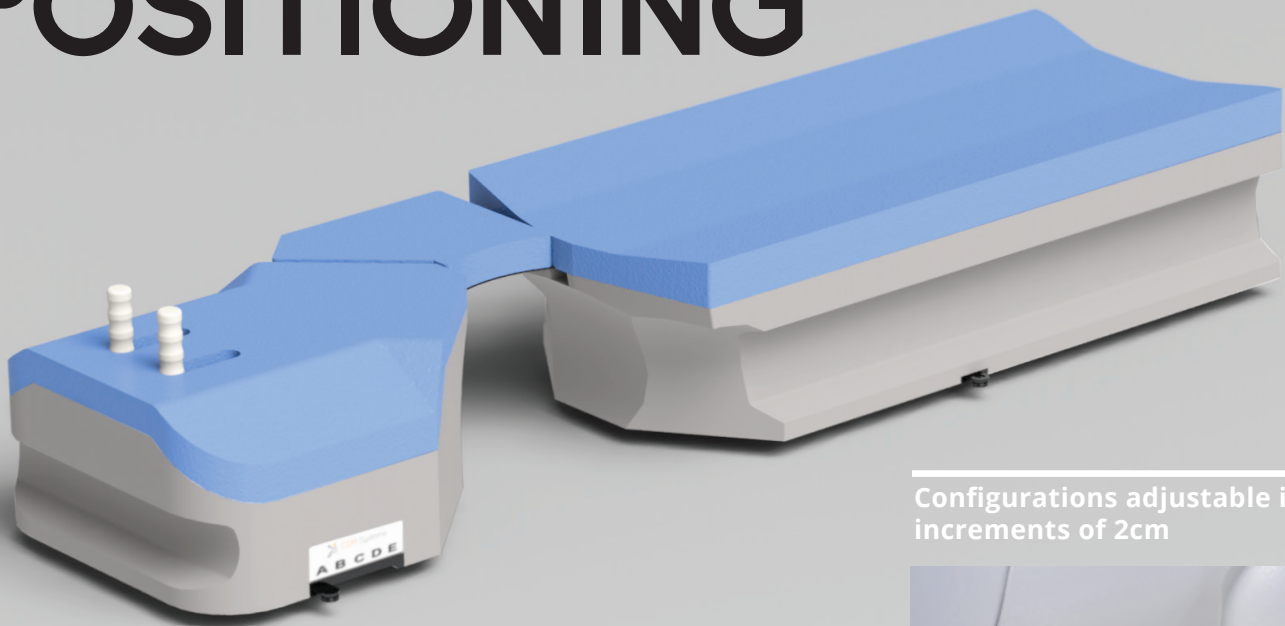

Configurations adjustable in increments of 2cm

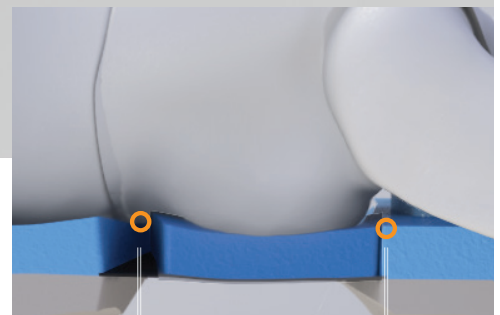

18cm configuration

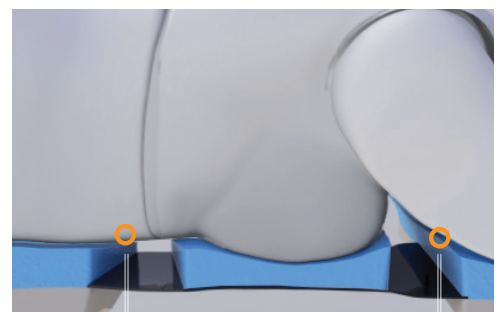


28cm configuration

PRODUCT OVERVIEW

Engineered for accurate quick setup of breast patients in the prone position. Allowing for lateral offsetting of patient when necessary.

The benefits of the Prone Breast technique have been well documented and include:


- Displacement of the lumpectomy cavity away from chest wall, heart and lung
- Reduction in cardiac systole and respiratory movement
- Enhanced cosmetic results with reduced soreness and skin reaction
- Optimized daily reproducibility

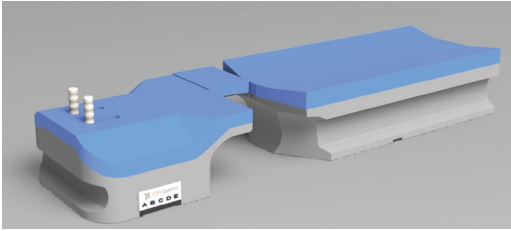

PB06

Standard height with standard cushion set Advanced Prone Breast Patient Positioning elevates chest wall 24cm from couch top, this design provides maximum patient support. The treatment breast opening is adjustable from 18-28cm.


PB06-2cm: Order ID to replace standard cushion set with firm cushion set.

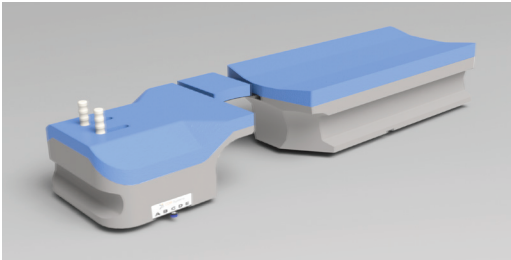

PB06-LP

Low profile Advanced Prone Breast Patient Positioning elevates chest wall 17 cm from couch top, this design provides maximum patient support. The treatment breast opening is adjustable from 18-28cm.


PB06-LP-2cm: Order ID to replace standard cushion set with firm cushion set.

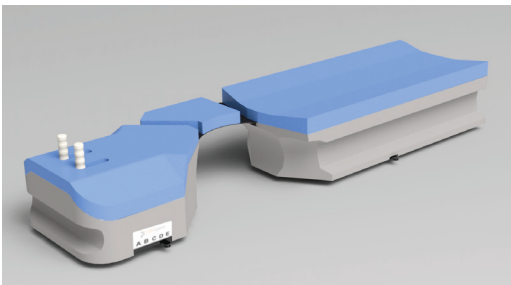

PB06-S

Standard height Advanced Prone Breast Patient Positioning elevates chest wall 24 cm from couch top. The open design is engineered for Superclav access and multi-field treatments. The treatment breast opening adjustable from 18-28cm.

PB06-S-2cm: Order ID to replace standard cushion set with firm cushion set.

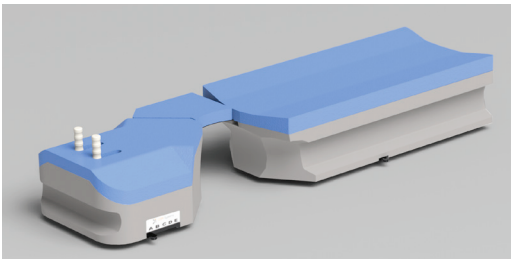

PB06-S-LP

Low profile Advanced Prone Breast Patient Positioning elevates chest wall 17 cm from couch top. The open design is engineered for Superclav access and multi-field treatments. The treatment breast opening adjustable from 18-28cm.

PB06-S-LP-2cm: Order ID to replace standard cushion set with firm cushion set.


MR For a MR Safe version of any of the Prone Breast boards, add a "K" to the order ID.

The Contralateral Breast Support Plate and cushion are Interchangeable for left or right side


Low Profile elevates patient 17cm form couch top for accommodating patients on all CT/ Sim including smaller bore scanners.

Standard height elevates patient 24cm from couch top specifically optimized for large bore CT/Sim.

Indexing bars sold separately.

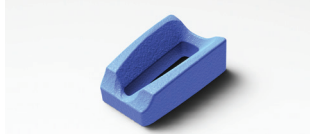

PF01
Prone Head Support Module

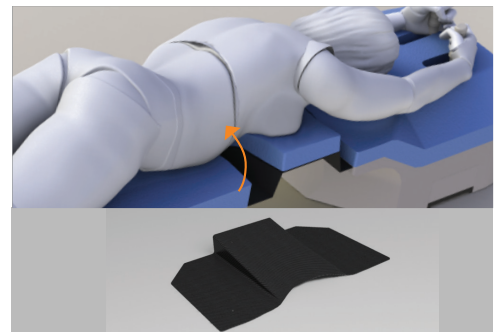

PB06-WDG

MR 10° angled Contralateral Breast support wedge compatible with all PB06 models and is adjustable for all opening setting of 18-28cm